

Portable Decontamination and Sterilization System

William Bell, Suzanne Smerjac and Bryan Smith

**Joint Service Scientific Conference on Chemical and
Biological Defense Research**

Hunt Valley, MD, 17 November 2004

TDA Research Inc. • Wheat Ridge, CO 80033 • www.tda.com

Report Documentation Page				Form Approved OMB No. 0704-0188	
Public reporting burden for the collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to a penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number.					
1. REPORT DATE 17 NOV 2004		2. REPORT TYPE N/A		3. DATES COVERED -	
4. TITLE AND SUBTITLE Portable Decontamination and Sterilization System				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S)				5d. PROJECT NUMBER	
				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) TDA Research Inc. Wheat Ridge, CO 80033				8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)				10. SPONSOR/MONITOR'S ACRONYM(S)	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release, distribution unlimited					
13. SUPPLEMENTARY NOTES See also ADM001849, 2004 Scientific Conference on Chemical and Biological Defense Research. Held in Hunt Valley, Maryland on 15-17 November 2004 . , The original document contains color images.					
14. ABSTRACT					
15. SUBJECT TERMS					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT UU	18. NUMBER OF PAGES 22	19a. NAME OF RESPONSIBLE PERSON
a. REPORT unclassified	b. ABSTRACT unclassified	c. THIS PAGE unclassified			

Abstract

TDA Research, Inc., (TDA) is developing a portable system to generate chlorine dioxide, which can be used for biodecontamination of small items and to sterilize medical and dental instruments in austere environments. Prototype systems can effectively sterilize both biological indicators and medical instruments inoculated with bacterial spores. At room temperature we found D-values (time for 90% kill) consistent with a sterilization cycle time of less than 30 minutes. The system has multiple features to assure simple and reliable operation. It requires no external power and includes indicators to verify correct functioning.

Organization of Presentation

- **Related work in decontamination at TDA**
- **Related work in portable devices**
- **Portable sterilization/decon system**
 - Overall design
 - Gas-generating reactions
 - Packaging of reactive ingredients
 - Indicators
 - Performance and applications

Decontamination Research at TDA

- **Protection of Building Occupants**
 - Gas-phase decontaminants, focus on bacterial spores
- **Spray-On Coating**
 - Facilitates decontamination of personal equipment
- **Catalytic Coatings**
 - Detoxification of CW agents under ambient conditions
 - HD oxidation catalysts
 - VX hydrolysis catalysts

Portable Devices at TDA

No Batteries or External Power Source Required

- **Self-heating products**

- Heater for military rations
 - Solid, activate by adding water
 - Designed to heat food pouches like US Army MRE
 - Being evaluated in US and other countries
- Self-contained package to heat liquids and oils
 - Easily activated by breaking frangible seal; product does not contact heater

- **Mosquito trap**

- Survey disease-carrying insects, monitor control programs
- Self-contained source of CO₂ attractant

- **Portable sterilizer/decon unit**

Sterilizer/Decon System

- **Sterilizer for medical/dental instruments, biodecon of small items**
 - For military medics, humanitarian relief workers
 - Medics can take fewer units of each instrument, sterilize for re-use, so fewer medical supplies must be carried
 - Sterilant generated by chemical reaction
 - No external power required
- **Gas-phase sterilant generated by chemical reaction**
 - Sterilant is chlorine dioxide, ClO_2
 - Rapid, proven effective in earlier studies
 - Solutions of ClO_2 in water are effective against VX and HD on surfaces; no data for gaseous ClO_2 with CW agents
- **ClO_2 is effective against bacterial spores**
 - Used to decontaminate anthrax spores in buildings
 - Bacterial spores are severe challenge, used in U.S. FDA certification protocol for sterilizers
 - Sterilization in less than 30 minutes

Sterilizer Unit and Subsystems

- **Unit is lightweight plastic bag**
 - Folding seal with adhesive closure, gas-tight
 - Similar to currently used resterilization pouches
 - Cannot easily re-open
 - Peelable foot for easy access after sterilization cycle
 - Obvious if opened, indicates sterilizer has been used
- **Subsystems:**
 - Gas-generating reactions
 - Packaging of reactive ingredients
 - Sterilizer bag design (single use)
 - Indicators
 - Vent with ClO_2 scavenger
 - Overwrap

Generation of ClO_2

- **Widely used in water purification**
- **Cannot be shipped, must be prepared at point of use**
- **Common ClO_2 production methods:**
 - Electrochemical processes
 - Sodium chlorate reduction with SO_2 or organics
 - Acidification of sodium chlorite


Gas-generating Reactions in Portable System

- **Generate mixture of ClO_2 and CO_2**
 - React solid acid and/or anhydride (HA) with solid NaHCO_3 and solution of NaClO_2 in water:




- **Use of gas mixture:**
 - Increases safety: high concentrations of ClO_2 (>10%) can explode, but our low concentration is stable
 - Improves gas distribution: larger total volume of gas produced
 - Adding salt to the aqueous solution used to produce sterilant allows humidity control, which improves kill of bacterial spores

Packaging Reactive Ingredients for Long Shelf Life

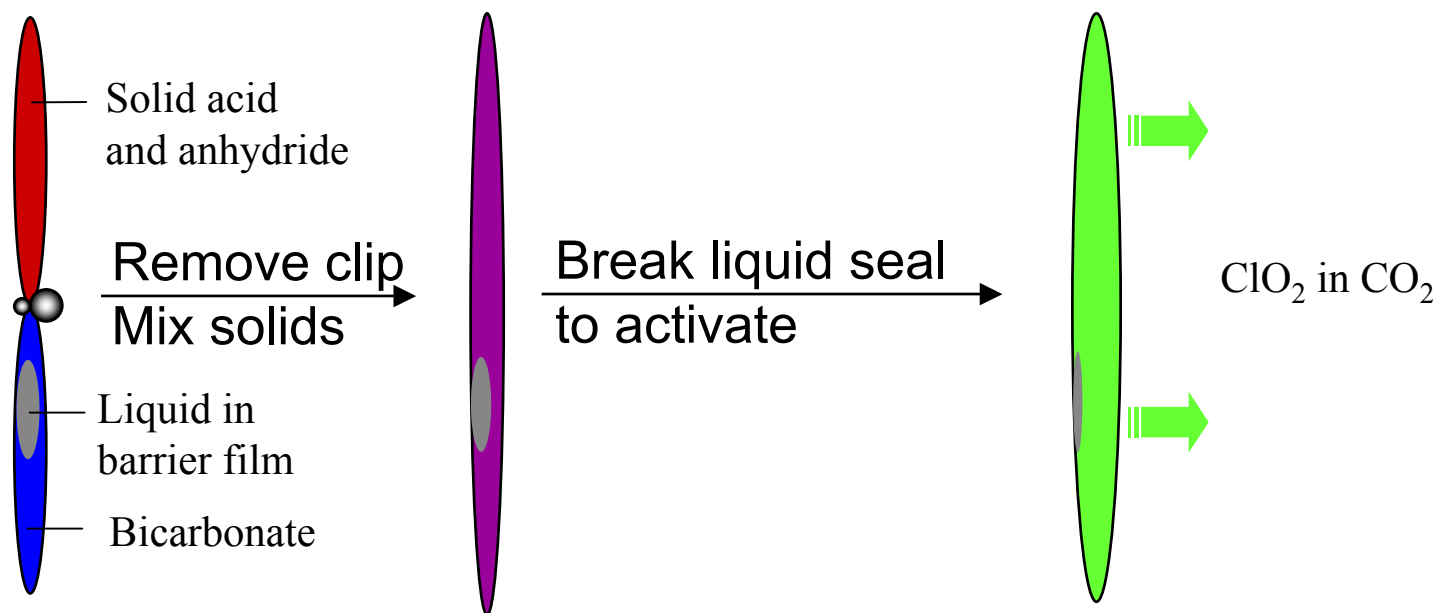


Schematic of Gas Generator

Side View

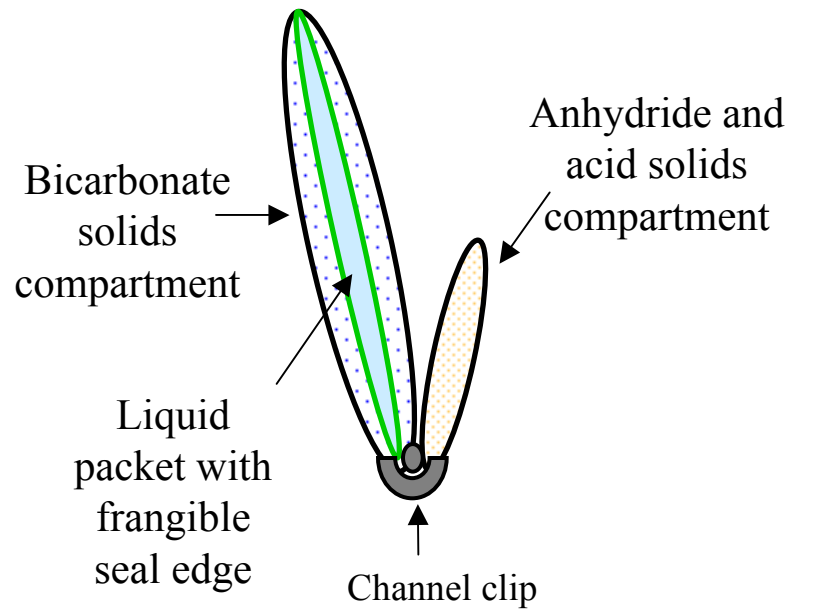


Liquid packet is attached to edge or side of chemical packet

Top view (unfolded)

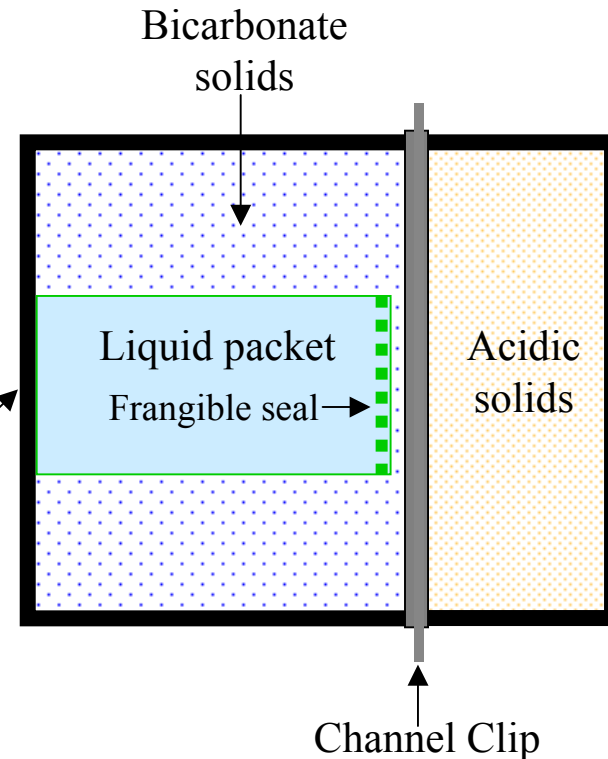

Gas Generator Performance


Chlorine Dioxide Concentration vs. Time

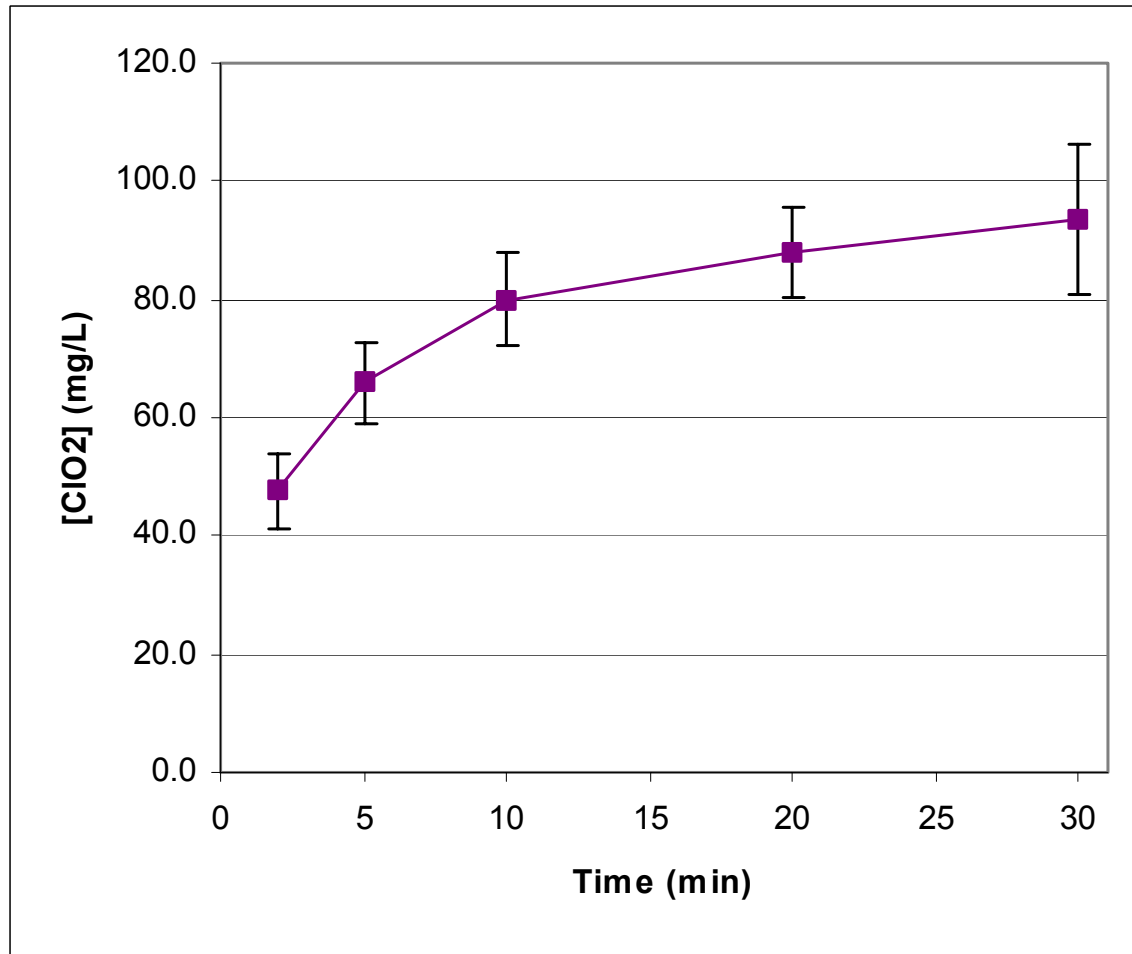

Sterilizer Materials (1)

- **Outer film**

- Transparent polyolefin, polyester, nylon, or laminates to see indicators inside
- Inner heat-seal surface
- Colored to avoid photodecomposition of ClO_2
- Best: polyester (dye uptake) laminated to polyolefin (heat sealing)

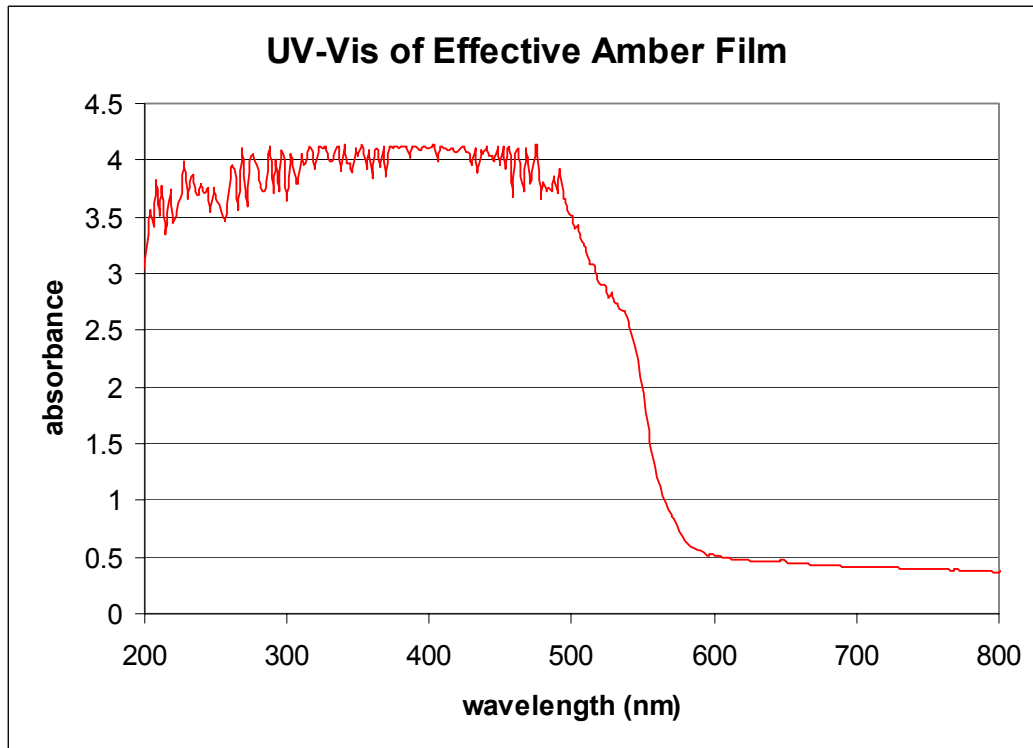


Film must block violet and UV light to prevent photodegradation of ClO_2

Indicators

- **Indicator to show ClO_2 is present**
 - Indicator is printed on or attached to inside of pouch
 - Changes color on exposure to ClO_2
 - Located at side of sterilizer; ClO_2 contacts this indicator after it has spread through the sterilizer
- **Temperature**
 - Reversible indicator on outside of sterilizer bag to show when unit is too cold and should be warmed for correct functioning

Vent with ClO_2 Sorbent

- **Use activated carbon to absorb ClO_2**
 - To assure operator safety
 - Carbon filter units simple, inexpensive
- **If filter is sufficiently far from sterilizer interior, no check valve needed**
 - Calculations indicate that 10 cm is far enough
- **Requires fitting to attach tubing to sterilizer bag and manually powered pump**

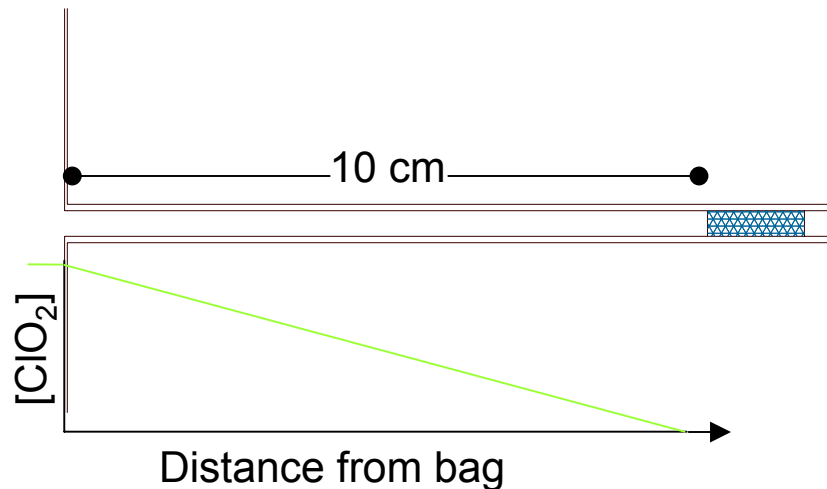


Test Methods

- **Use Biological Indicators (BIs)**
 - Each BI is 1-cm diameter stainless steel disc, inoculated with *ca.* 10^6 spores of *B. atrophaeus*, in Tyvek/Mylar bag.
- **Test in plastic bag; generate ClO_2 from solid mixture**
 - During sterilization, ClO_2 penetrates Tyvek
- **Vary exposure time, ClO_2 concentration, temperature, humidity**


D-Value Calculation

- D-value (time for 1-log population reduction) from population data (graph below) is 1.6 minutes
- Growth/no-growth analysis by Stumbo-Murphy-Cochran equation indicates D-value of 1.1 minutes
- Also effective when instruments are in standard Tyvek resterilization pouches

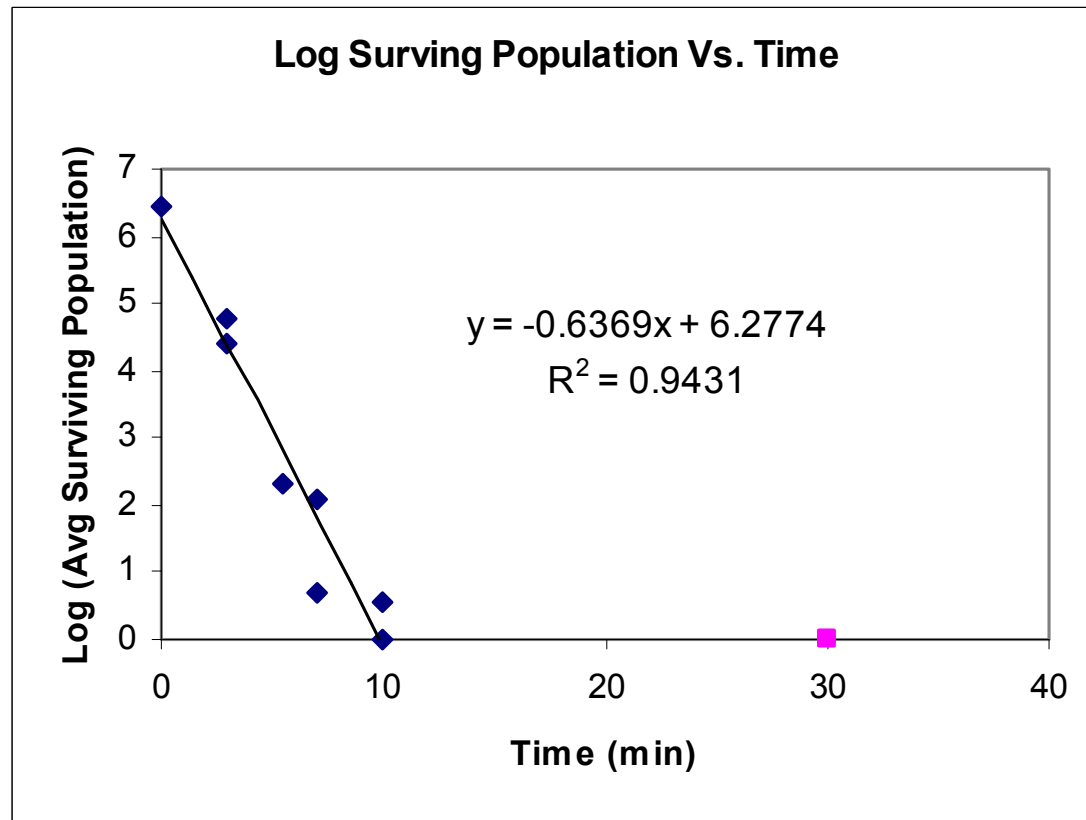

Humidity Control

To eliminate bacterial spores, it is highly desirable to control the humidity during sterilization to 70 to 95%, preferably 90 to 95%. Adding NaCl to the NaClO_2 solution lowers the humidity of the gases produced from 100% to the desired range.

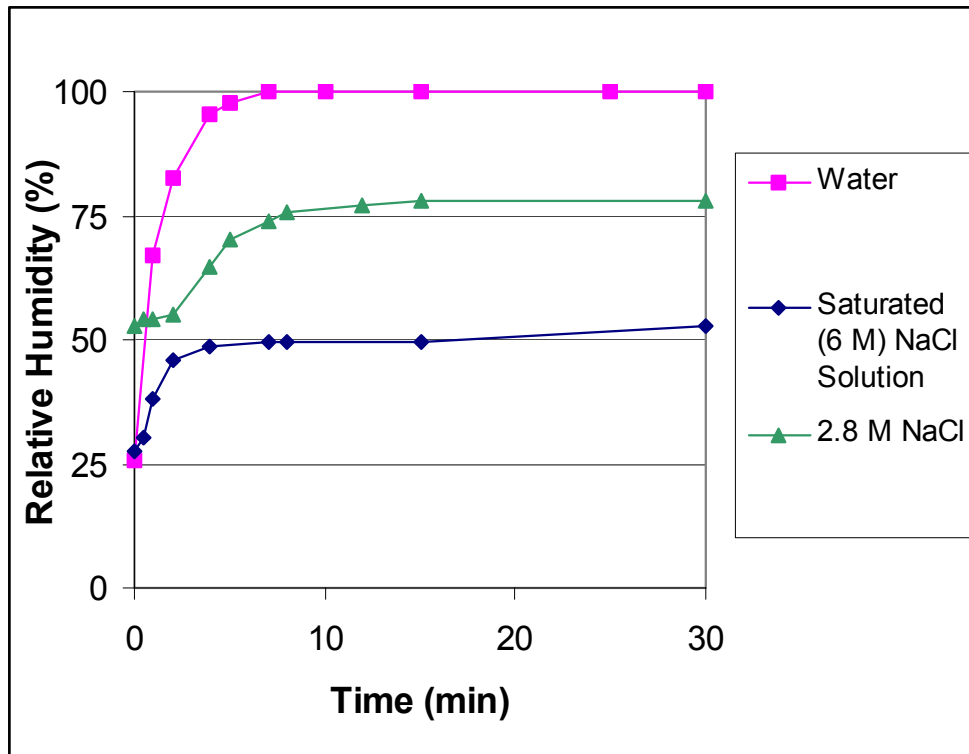

- In a single-use collapsible bag, most of the gas in the bag during sterilization is from the gas generator.
- By preventing condensation, we minimize the time for decontamination and help prevent corrosion

Performance and Dimensions

- **Demonstrated generation of effective ClO₂ dose and humidity control in portable device**
- **Tested with biological indicators (BIs), and with instruments inoculated with bacterial spores**
- **Sterilization in less than 30 minutes at room temperature**
 - Longer treatment times required at low temperature
 - Reduced efficacy when spores were deposited in combined serum/hard water matrix; instruments must be thoroughly cleaned before treatment to assure sterilization
- **Size (largest tested to date)**
 - 56 cm by 36 cm
 - Weight 125 g
- **Can process medical instruments weighing more than 10 times as much as the sterilizer**

Summary

- **Convenient source of ClO₂ gas**
- **Readily transported**
- **Requires no batteries or external power**
- **Suitable for sterilization of medical/dental instruments, biodecontamination of small items**
- **TDA has applied for patent**
- **Development is continuing**

Acknowledgement and Disclaimer

- **This work was supported by the U.S. Army Medical Research and Materiel Command under Contract No. DAMD17-03-C-0023.**
- **The views, opinions and/or finding contained in this report are those of the authors and should not be construed as an official Department of the Army position, policy or decision unless so designated by other documentation.**

References

Gates, D. (1998). The Chlorine Dioxide Handbook. Denver: American Water Works Association.

**Haas, C.N. (2001). “Decontamination Using Chlorine Dioxide.” Testimony for the Committee on Science, U.S. House of Representatives, Nov. 8.
www.house.gov/science/full/nov08/haas.htm.**

Harper, B., L. Larsen, and P. Stevenson (2002). “Laboratory Validation of Chlorine Dioxide Decontamination (Anthrax Spores),” presented at the Joint Service Chemical and Biological Decontamination Conference, San Diego, CA, 22-24 October.

Leighton, T., and G. Eggum and K. Wheeler (2002). “Use of Chlorine Dioxide for Large-Area Decontamination” presented at the Joint Service Chemical and Biological Decontamination Conference, San Diego, CA, 22-24 October.

Leighton, T, K. Wheeler, L. Wein, S. Barry and R. Brookmeyer (2004). “Decontamination of Bacterial Threat Agents with Chlorine Dioxide.” Presentation at DECON 2004 conference, Tampa, FL, May 17-20.